




Miele

EasyLoad – Fast loading and highly efficient use of chamber space.
The new loading system for cleaning laboratory glassware.

Miele Professional. Immer Besser.

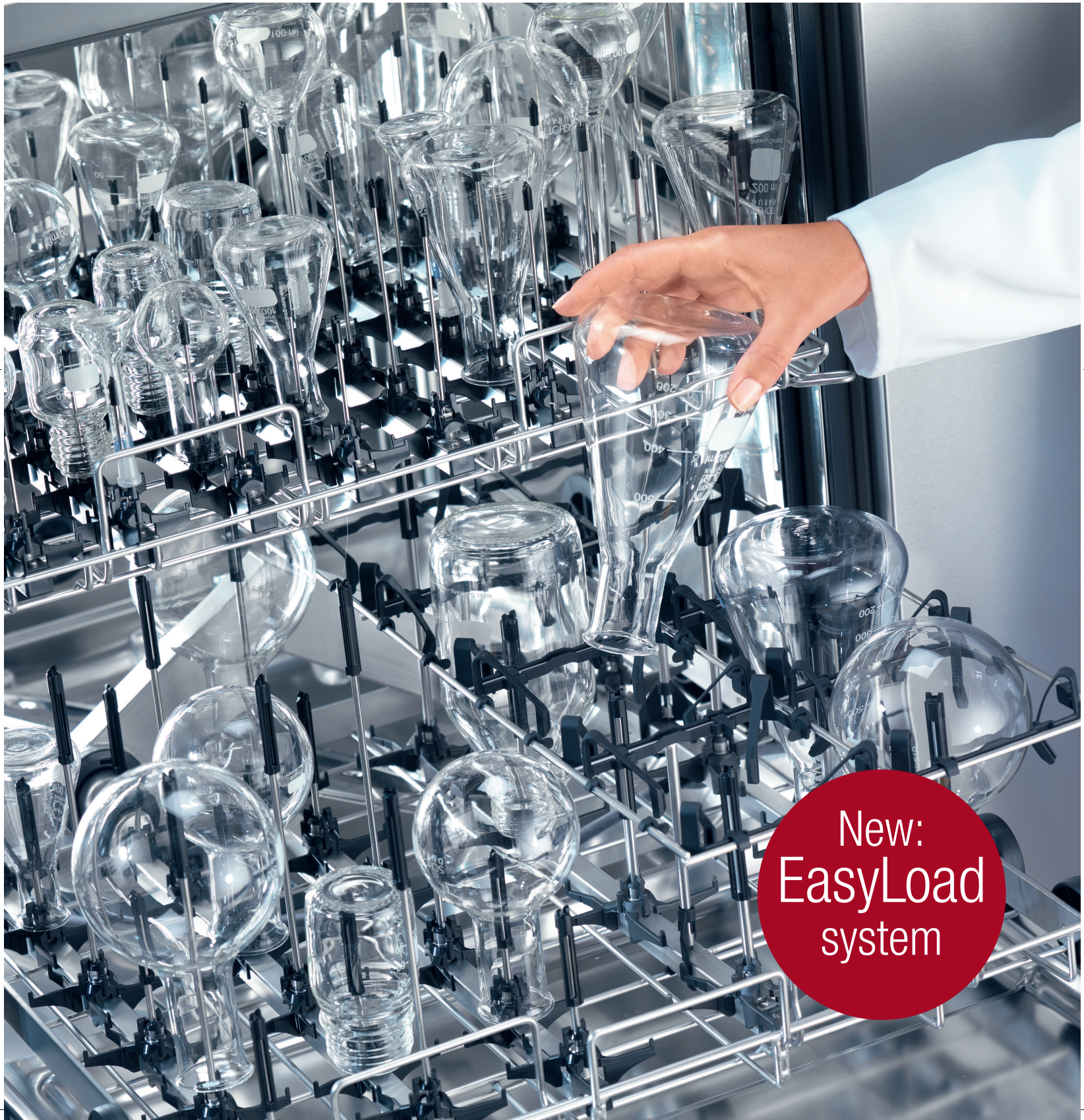

New:
EasyLoad
system

Miele reinvents loading: with EasyLoad.

- New injector system for all common injector modules and large number of further load carriers
- Optimum use of available vertical clearance
- Particularly fast and simple arrangement of load
- Thorough internal cleaning of glassware of all shapes and sizes
- Secure hold and exceptional material protection during wash cycle

Figure: Load rack with support frame




EasyLoad means efficiency and convenience across the board:

- New injector nozzles (3 diameters, 6 lengths)
- New plastic holders and rests
- New injector modules (8–32 injectors)
- New support frame (for 8–32 injector positions)
- Available for PG 8583 and PG 8583 CD, PG 8593



The use of the new injector nozzles and rests without a support frame is possible in virtually all Miele lab load carriers.

Even if a nozzle tip rests against the base of a glass, the new injectors still guarantee thorough internal cleaning. The new rest surrounds the injectors in the module or basket, preventing lateral movement without wasting valuable height.

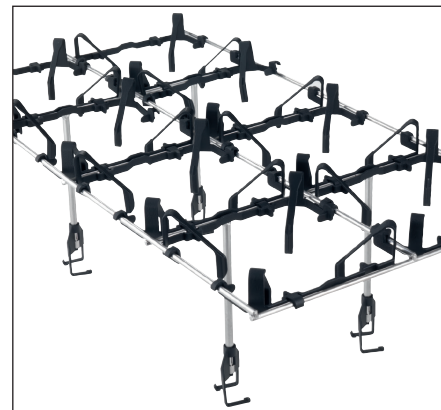

The new support frame provides additional hold. As and when required, it is placed on top of the load carrier, providing a secure hold to even large items of laboratory glassware in a vertical and protected position.

EASYLOAD



System4Lab